



TECNILABO

CELLCARE

CO₂ INCUBATOR
EXCELLENCE FOR YOUR LAB



AVAILABLE BY 2024


PRECISION & COMFORT

- The Heart - ensure stable humidity and CO₂ levels and accurate readings with a multi-parameter sensor
- O₂ control and regulation system available on request
- The Brain - keep your cells safe with an advanced electronic control system
- Real-time data storage with a user-friendly touch screen interface with USB and ethernet connectivity
- No internal spare parts to handle (e.g. UV lamps or HEPA filter)



CLEANABILITY & CONTAMINATION CONTROL

- 180° C (356° F) decontamination
- Ensure optimal cell growth thanks to a seamless deep drawn chamber with fanless design
- Removable and autoclave-safe shelving system and water tray
- Smooth external surfaces with painted steel

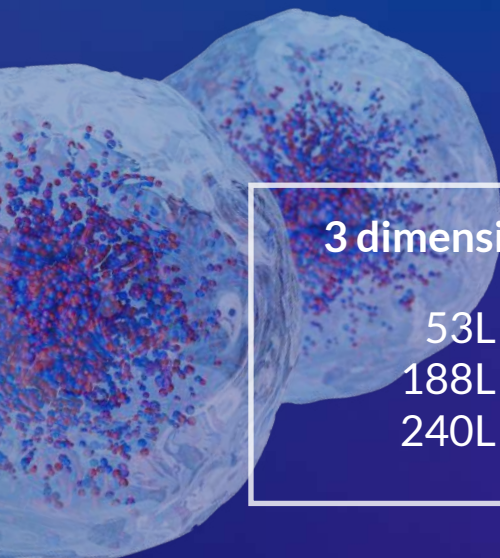


3 dimensions

53L
188L
240L

PERFORMANCE & SPEED

- Proprietary advanced direct heating system for temperature uniformity
- CO₂ and temperature recovery <5 minutes
- Protect cells from vibration and turbulence with a seamless design

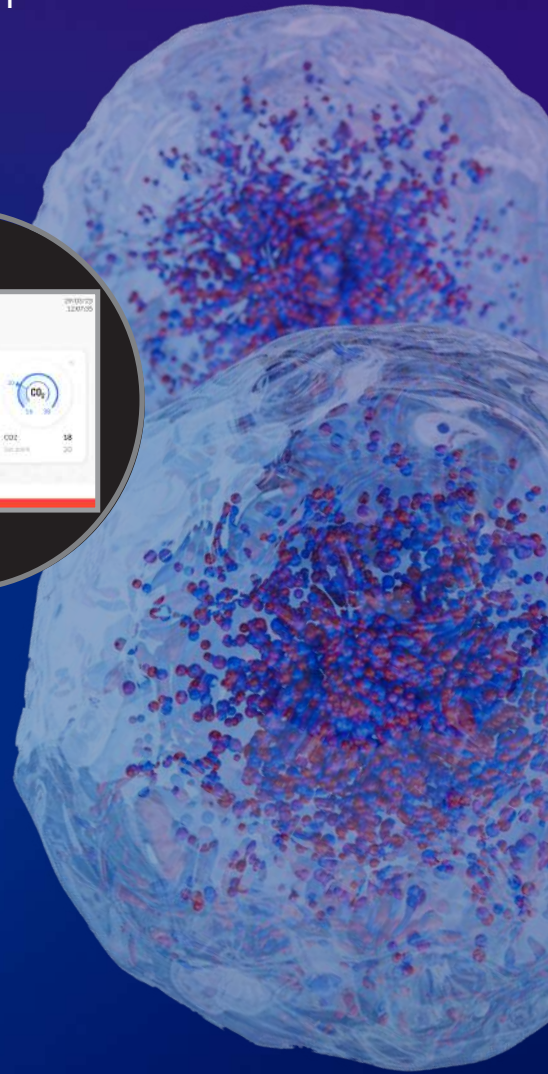

SAFETY

- No need to remove sensor during high temperature cycle
- Anti-vibration soft door with available automatic locking system

DIMENSIONS

	MAXI	MIDI	MINI
External dimensions (in)	Coming 2025	27.79(W) x 30.31(D) x 43.9 (H)	24.01(W) x 24.41(D) x 37.07(H)
Internal dimensions (in)	Coming 2025	20.87(W) x 19.68(D) x 27.16 (H)	15.75(W) x 13.78 (D) x 15.75 (H)
Volume (L)	240	188.6	53

FIND OUT MORE



YOUR SAFETY
IS OUR COMMITMENT



TECNILABO HEADQUARTERS

TECNIPLAST S.P.A.
VIA I MAGGIO, 6
21020 BUGUGGIATE (VA)
tel +39 0332 809711 - fax +39 0332 458315
info@tecnilabo.com

SISTER COMPANIES

AUSTRALIA/NEW ZEALAND | Tel. + 61 2 8845 6500 • www.tecniplast.it • info@tecniplast.com.au
CANADA | Toll Free: +1 855.347.8718 • www.tecniplastcanada.com • info@tecniplastcanada.com
CHINA | Tel. + 86 (0)21 50810920 • www.tecniplast.cn • tecniplastchina@tecniplast.it
FRANCE | Tél. +33 (0)4 72 52 94 41 • www.tecniplast.fr • info@tecniplast.fr
SWITZERLAND | Tel +41 (0)81 750 31 40 • www.tecniplast.ch • info@tecniplast.ch
UNITED KINGDOM | Tel. + 44 (0)345 0504556 • www.tecniplastuk.com • info@tecniplastuk.com
USA | Toll Free: +1 877.669.2243 • www.tecniplastusa.com • info@tecniplastusa.com